

Intensiv Ortho-Strips® System



Oscillating diamond strips for the double-sided or one-sided interdental reduction of the tooth size in orthodontics

When opening the interdental space and reducing, finishing and polishing enamel in orthodontics, care must be taken to avoid grooves and scratches. It is also important to proceed systematically, following a coarse-to-fine sequence in selecting grits. Unlike manual strips, the Intensiv Ortho-Strips System provides rapid, controlled enamel reduction followed by polishing without unnecessary removal of healthy tooth substance.



Product Description

- Flexible strips, double-sided or one-sided diamond coating
- 6 different grits:
 - 08 µm, orange, opening
 - 40 µm, red, contoured reduction
 - 60 µm, brown, fast reduction
 - 80 µm, blue, very fast reduction
 - 25 µm, white, finishing
 - 15 µm, yellow, polishing
- Sterilizable

Intensiv Ortho-Strips System, Double-Sided, assortment of 6 strips, REF OSSet02/6

- OS08OP-DS, OS40M-DS, OS60C-DS, OS80XC-DS, OS25F-DS, OS15POL-DS

Intensiv Ortho-Strips System, Left-Sided, assortment of 6 strips, REF OSSet02L/6

- OS08OP-L, OS40M-L, OS60C-L, OS80XC-L, OS25F-L, OS15POL-L

Intensiv Ortho-Strips System, Right-Sided, assortment of 6 strips, REF OSSet02R/6

- OS08OP-R, OS40M-R, OS60C-R, OS80XC-R, OS25F-R, OS15POL-R

Intensiv Ortho-Strips System, Opener

- Strip with serrated edge and ultra-fine 8µm grit to open contact points without friction



Indications

- Opening of the interdental contact points
- Widening of the interdental spaces in orthodontics by bilateral or unilateral tooth size reduction
- Elimination of minor crowding and treatment finishing in orthodontics
- Bilateral proximal enamel polishing

Benefits

- Efficient opening of the interdental contact points
- Rapid and controlled enamel reduction
- Proximal contouring, finishing and polishing of both adjacent teeth in a single procedure
- Reusable

Dr. Francesco Garino, Turin, Italy



Clinical pictures



- 01 Opening of the interdental space with OS08OP-DS, Opener, 8 µm
- 02 Contoured reduction with OS40M-DS, Medium, 40 µm
- 03 Fast reduction with OS60C-DS, Coarse, 60 µm
- 04 Very fast reduction of enamel with OS80XC-DS, Extracoarse, 80 µm
- 05 Controlled finishing with OS25F-DS, Fine, 25 µm
- 06 Polishing with OS15POL-DS, Polishing, 15 µm

Colour code	Micron	Description
■ 524	80 µm	Extracoarse
■ 524	60 µm	Coarse
■ 514	40 µm	Medium
■ 514	25 µm	Fine
■ 504	15 µm	Polishing
■ 494	8 µm	Opener

Diamond-coated:
DS = Double-Sided
L = Left
R = Right

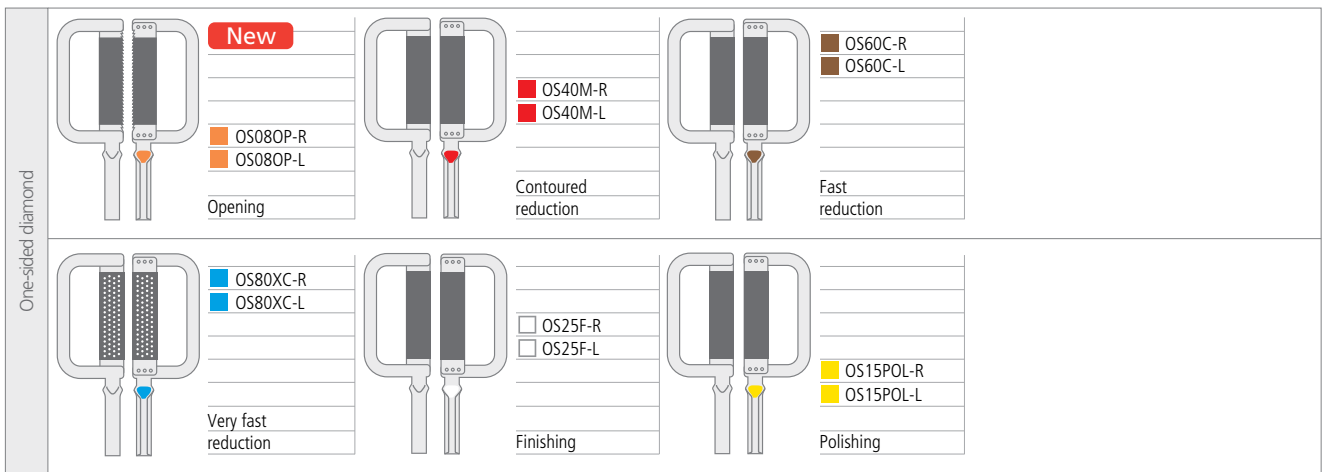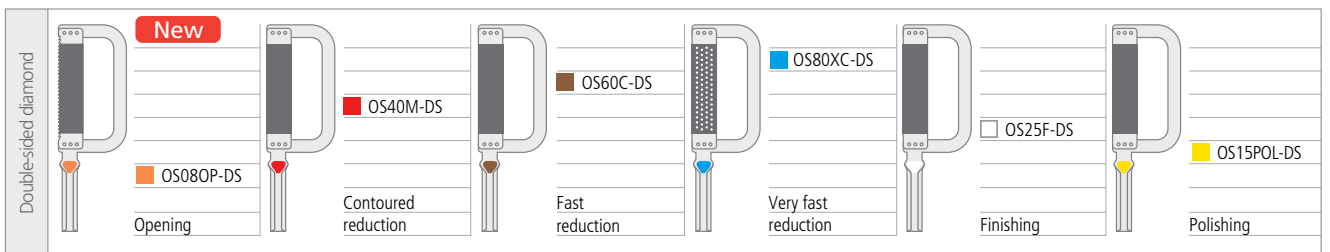

Intensiv Ortho-Strips Double-Sided and One-Sided are available in packages of 3, 6 or 12 pcs.