

Intensiv Swingle[®] Pat. EP 2754406B1

SWISS EDITION

Reciprocating contra-angle handpiece with micromechanical dual function (axially free-rotating / axially fixed position of the oscillating instruments) - Durable, refined and efficient, including additional service package

Efficient use of a reciprocating contra-angle handpiece requires a choice between the axially fixed position and the free axial rotation of the oscillating instruments during the same treatment. Intensiv Swingle SWISS EDITION with a better grip design and a scratch and wear-resistant glass-coating. This special coating provides the contra-angle handpiece with long lasting protection against scratched surfaces.



REF WG-69LTSWISS
(with light)

REF WG-69ASWISS
(without light)



Toggle switch fixed downwards: oscillating instruments in axially fixed position
Toggle switch fixed upwards: axially free-rotating oscillating instruments.

Product Description


DUAL INTENSIV SWINGLE CONTRA-ANGLE

- Better grip design of the contra-angle handpiece
- Special coating provides lasting protection against scratched surfaces
- Reciprocating contra-angle handpiece with integrated water spray and optional light
- Stroke movement 0.9 mm
- 20,000 oscillating stroke movements at 40,000 rpm
- With toggle switch to select axially free-rotating or fixed position of the oscillating instruments

A toggle switch on the contra-angle head allows you to choose between axially free-rotating or axially fixed oscillating instruments.

Indications

- Swingle: for oscillating instruments with oscillating movement
- Swingle/Proxoshape: Precise preparation and margin finishing
- Swingle/Proxoshape: Fast and efficient excess removal, finishing and polishing of restorations
- Swingle/Ortho-Strips: Approximal enamel reduction, IPR, in orthodontics

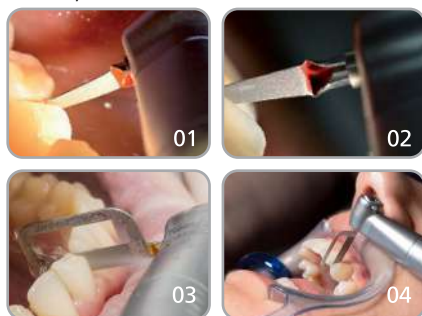

01-02) Dr. A. Devigus, Bülach, Switzerland
03-04) Dr. F. Garino, Turin, Italy

Benefits

- Lasting integrity of the contra-angle's surface, thanks to the scratch resistant and durable special coating
- 20.000 oscillations per minute for better use performance of files
- 3 year Warranty including one maintenance service

Clinical pictures

- 01 Approximal finishing with Intensiv Swingle SWISS EDITION and Intensiv Proxoshape Flexible, red, 40 µm
- 02 Finishing of margins with Intensiv Swingle SWISS EDITION and Intensiv Proxoshape Flexible, red, 40 µm
- 03 Approximal contact point opening with Intensiv Swingle SWISS EDITION and Ortho-Strips Opener, orange, 8 µm
- 04 Approximal enamel reduction with Intensiv Swingle SWISS EDITION and Ortho-Strips Central, red, 40 µm